


진단이



1
Sampling

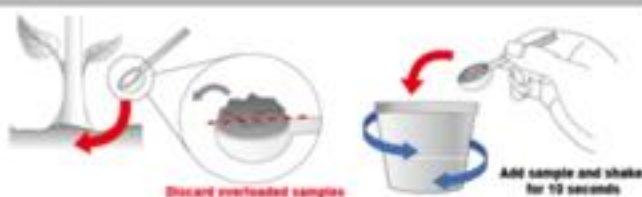

2
Use the Supernatant layer



3
Loading



4
Reading



5
Discard



Interpretation	Line Intensity	(C vs T)	Fusarium Titer	Ralstonia Titer
Low		C>>>T	$\leq 10^3$ cfu/ml	$\leq 10^4$ cfu/ml
Low (Early stage)		C>T	$\leq 10^4$ cfu/ml	$\leq 10^5$ cfu/ml
Middle (Middle stage)		C=T	10^5 cfu/ml	10^6 cfu/ml
High (Late stage)		C<<<T	$\geq 10^6$ cfu/ml	$\geq 10^7$ cfu/ml

DISTRIBUIDOR EXCLUSIVO

KIT RAPIDO PARA DETERMINACION DE ENFERMEDADES DEL SUELO Y LA PLANTA

MADE IN



KOREA



SUMIKIT

DETERMINACION

- FUSARIUM OXYSPOURUM
- RALSTONIA SOLANACEARUM
- ERWINIA
- PHYTOPHTORA
- ✓ • NEMATODOS
- ✓ • CMV VIRUS
- ✓ • Cgmmv
- ✓ • TSWV

NEW

2

10 Min resultado in sitio

3

99.9% DE EFICACIA

Gracias a su tecnologia patentada, es seguro y de facil uso

4

Ensayo Inmunocromatografico en fase solida para la detección rápida y semicuantitativa de hongos en extractos vegetales.